




Case-by-Case
MOBILE SERVICE

Storz Modulith SLX-F2

Portable Lithotripsy System

Ideal for custom treatment strategies in which high pressure and energy concentration are necessary.

The Storz Modulith® SLX-F2 is designed for the treatment of renal and ureteral calculi using extracorporeal shock wave lithotripsy (ESWL).

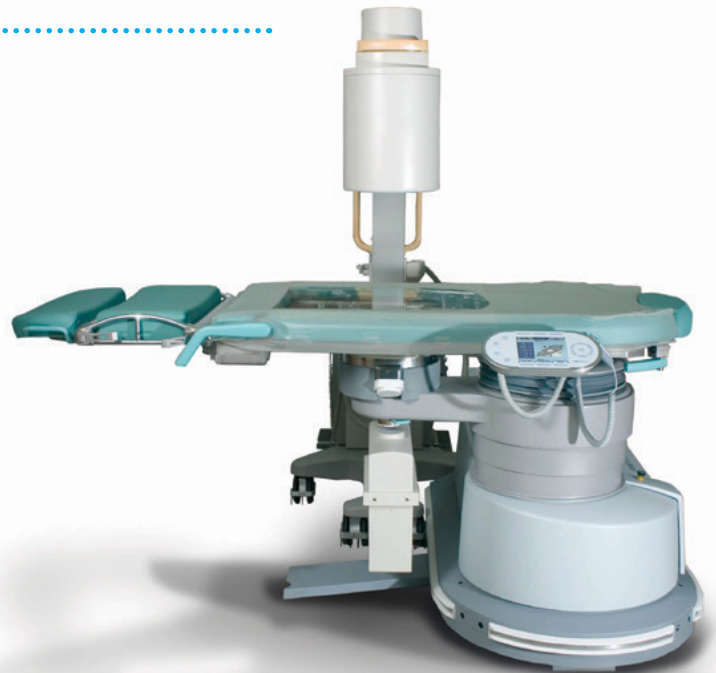

[THE FORTEC ADVANTAGE]

- ✓ Avoid capital and maintenance costs
- ✓ Experienced, professional technician support
- ✓ A stable, professional partner
- ✓ Large inventory of surgical technologies

Call today

800.963.7101

to start ForTec Mobile Service

Storz Modulith SLX-F2

Portable Lithotripsy System

Using the Modulith SLX-F2, it is possible to have a customizable treatment option based on your patient needs and anatomical conditions. High pressure and precise energy concentration are often necessary in the disintegration of ureteral stones, whereas lower pressures and a wider focal zone can be effectively used for larger stones located in the kidney.

Versatility for Custom Treatment

Since no two stones or anatomies are the same, the SLX-F2 dual focus design provides a highly adaptable shock wave parameter for a completely customizable treatment. The dual focus allows you to customize a treatment strategy that best suits your patient's needs.

Low Side Effects

The effective and efficient cylindrical shock source of the Modulith SLX-F2 provides unsurpassed fragmentation efficiency with low side effects and a focal depth of 165 mm.

Technical Specifications:

Shock wave principle:	Electromagnetic dual focus cylindrical source
Penetration depth:	0-165 mm
Shock wave focus:	Dual focus 6 x 28 mm or 9 x 50 mm (size selectable during treatment)
Lithotripter dimensions:	169 x 89 x 124 cm
Lithotripter weight:	659 kg
Localized projections:	In-line localization
Ultrasound system:	Black and white (color Doppler unit)